
nedis

MANE110GN

Relaxing Massage Pillow

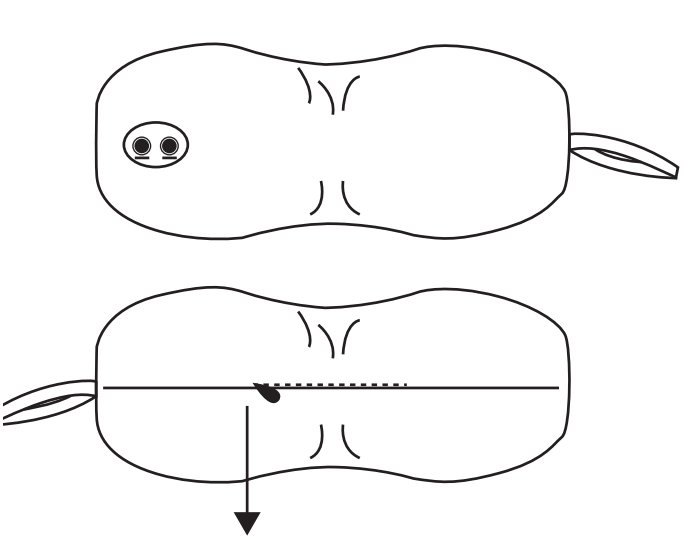Description

This product can effectively massage muscles, release arteries and veins and increase the blood circulation. It relieves stiffness, tiredness and paralysis while it can be used wherever: at home or at the office.

Product Features

- 6 massage modes (Ref. to table 1).
- Automatically switches off after 20 minutes.
- Improve the blood circulation of your head and neck for a better oxygen supply to the brain.
- The massaging device can be installed inside a pillow which can adapt to the different postures of the users and provide good support. This supportive function minimizes the harm caused by pressure and protects the backbone.
- This pillow is designed to fit a human's head curve ergonomically.

Mode	Massage Type
MODE1	Continuous gentle massage.
MODE2	Continuous strong massage.
MODE3	Changes from strong massage to stop massaging (Changes 70 times within 1 minute).
MODE4	Alternates gentle and strong massages. Changes between the two modes 60 times in 1 minute.
MODE5	Circulate massage mode 1, 2, 3 and 4 for around 3 seconds each.
MODE6	Circulate massage mode 1, 2, 3, 4, and 5 for around 45 seconds each.

Operation

- Press the ON/OFF button to switch on the massage pillow and activate the massage mode.
- Press the MODE button to select the different massaging modes.

Battery Replacement

When the massage pillow's vibration is slower or stops completely, the battery is low. Please replace the battery with 2pcs of AA batteries.

- Unzip the back of the pillow and take out the battery compartment.
- Press on the battery cover slightly and slide it down according to the arrow's direction on the cover.
- Replace the old batteries by the new ones.
- Close the battery cover. Put it back inside the massage pillow and zip the massage pillow again.

Attention

- Please put the massage pillow away from children.
- Do not use or store the massage pillow near magnetic fields and keep it away from high temperatures and moisture.
- Please take out the batteries and store the massage pillow properly when it is not going to be used for a long period.
- People without the ability to express themselves, such as babies, are forbidden to use this product.
- Do not throw the batteries into a fire.

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