
KN-BS10

Baby scale

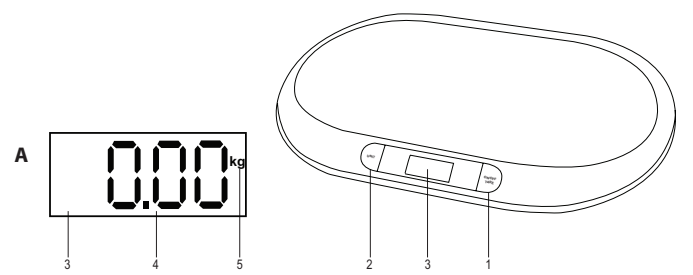Description

The baby scale is used to measure the body weight of babies and toddlers. (max. 20 kg).

Device (fig. A1)

1. On/Off button Tare button	- Press the button to switch on or off the device. - If you use a towel to lay the baby on, put the towel on the platform and press the button to set the display to "0.00".
2. Unit-of-measurement button	Press the button to adjust the unit of measurement: kilogram (kg), pound (lb) or stone (st).
3. Display	
4. Weight indicator	The indicator shows the body weight.
5. KG/LB/ST indicator	The corresponding indicator comes on when kilogram (kg), pound (lb) or stone (st) is set as unit of measurement.

Batteries

- Open the battery compartment.
- Insert the batteries (2x AA) into the battery compartment.
- Close the battery compartment.

Use

- Place the device on a stable, flat surface.
- To switch on the device, press the on/off button.
- Gently lay the baby or the toddler on the platform. The display shows the body weight.
- Remove the baby or the toddler from the platform.
- To switch off the device, press the on/off button.
- Note: The device switches off automatically after 20 seconds.

Safety

- Read the manual carefully before use. Keep the manual for future reference.
- Only use the device for its intended purposes. Do not use the device for other purposes than described in the manual.
- Do not use the device if any part is damaged or defective. If the device is damaged or defective, replace the device immediately.

Battery safety

- Use only the batteries mentioned in the manual.
- Do not install batteries in reverse polarity.
- Do not short-circuit or disassemble the batteries.
- Do not expose the batteries to water.
- Do not expose the batteries to fire or excessive heat.
- Batteries are prone to leakage when fully discharged. To avoid damage to the product, remove the batteries when leaving the product unattended for longer periods of time.
- If liquid from the batteries comes into contact with skin or clothing, immediately rinse with fresh water.

Cleaning and maintenance

Warnings!

- Do not use cleaning solvents or abrasives.
- Do not clean the inside of the device.
- Do not attempt to repair the device. If the device does not operate correctly, replace it with a new device.

Clean the outside of the device using a soft, damp cloth.

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