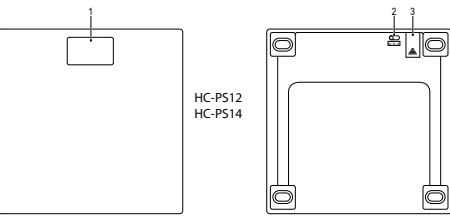
# KONIG

## HC-PS12 / HC-PS14

## Body scale

<span></span>	<span></span>	<span></span>
<b>A1</b>	HC-PS12	HC-PS14
		

## Description

## Device

1. Display	
2. Unit-of-measurement button	• Press the button to adjust the unit of measurement: kilogram (kg), pound (lb)
3. Battery compartment	• Open the battery compartment. <ul style="list-style-type: none"><li>HC-PS12 / HC-PS14</li></ul> <ul style="list-style-type: none"><li>Insert the battery (CR2032) into the battery compartment.</li></ul>

## Display indicators

"Lo"	Low battery	Replace the batteries.
"Er"	Overload	The person on the platform exceeds the maximum capacity of the body scale. Step off the platform to prevent damage.
"C"	Error	Step off the platform. Step on the platform to weigh yourself again.

## Use

- Place the device on a stable, flat surface.
- Switch on the device, press the platform with your foot.
- Step on the platform. Place your feet evenly on the platform without moving. The display shows your weight.
- Step off the platform. The device switches off automatically.

## Safety

- Read the manual carefully before use. Keep the manual for future reference.
- Only use the device for its intended purpose. Do not use the device for other purposes than described in the manual.
- Do not use the device if any part is damaged or defective. If the device is damaged or defective, replace the device immediately.

### Battery safety

- Use only the batteries mentioned in the manual.
- Do not install batteries in reverse polarity.
- Do not short circuit or disassemble the batteries.
- Do not expose the batteries to water.
- Do not expose the batteries to fire or excessive heat.
- Batteries are prone to leakage when fully discharged. To avoid damage to the product, remove the batteries when leaving the product unattended for longer periods of time.
- If liquid from the batteries comes into contact with skin or clothing, immediately rinse with fresh water.

### Cleaning and maintenance

- Do not use cleaning solvents or abrasives.
- Do not clean the inside of the device.
- Do not attempt to repair the device. If the device does not operate correctly, replace it with a new device.

- Clean the outside of the device using a soft, damp cloth.

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